

Using the Log Viewer


You can use the Log Viewer to look at the pressure trace of all the RAW files you acquired. All that data is stored on the on-board hard drive on the EASYnLC.

Using the Log Viewer	1
Copy files from EASYnLC hard drive to local PC	1
Compress files into a .zip file	2
Open .zip file with LogViewer.....	3
Mount the EASYnLC hard drive to file explorer:	5


Copy files from EASYnLC hard drive to local PC

Open file explorer and click on “admin (\\172.16.0.104) (Z:)”, if it prompts you for a log in and password use “admin” for both.


Note if the “admin (\\172.16.0.104) (Z:)” is not present or if you have to map for the first time follow the instructions at the end of this document.

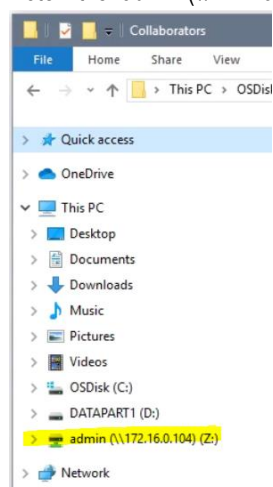

Double click on “batchLog”:

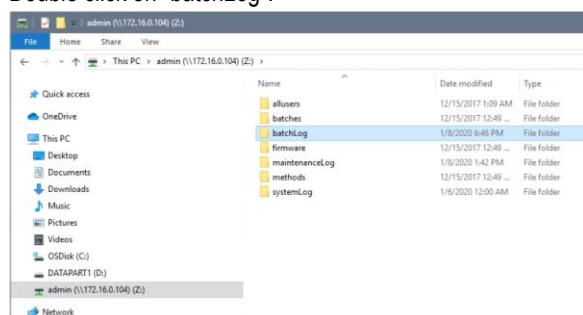

Now you should see a folder for each RAW file you acquired. Sort by "Date modified" to list the newest ones on top.

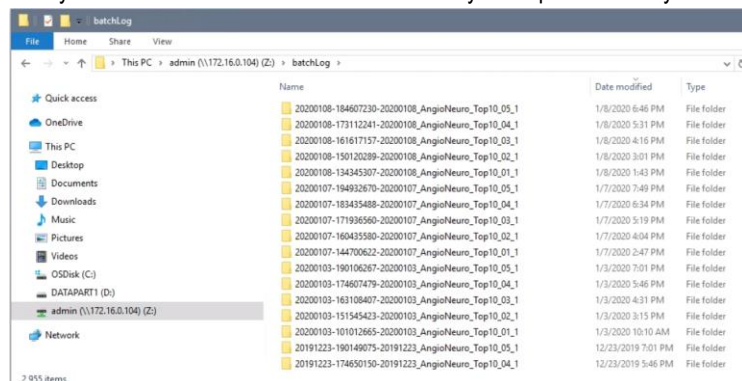


To copy the folders, drag the ones you want to look at to your folder on the local PC.

If you get a security error, acknowledge it by clicking "OK":




The folder should copy to the folder you selected:

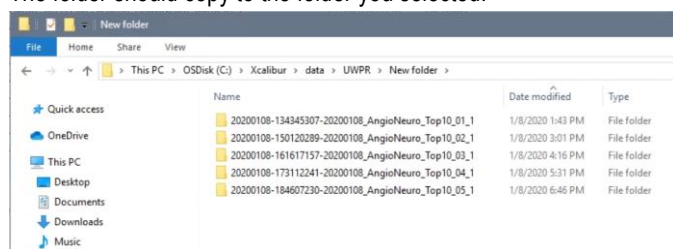


Compress files into a .zip file

To create a zip file, select all the folders, do a right mouse click and select "Send to" and then "Compressed (zipped) folder":

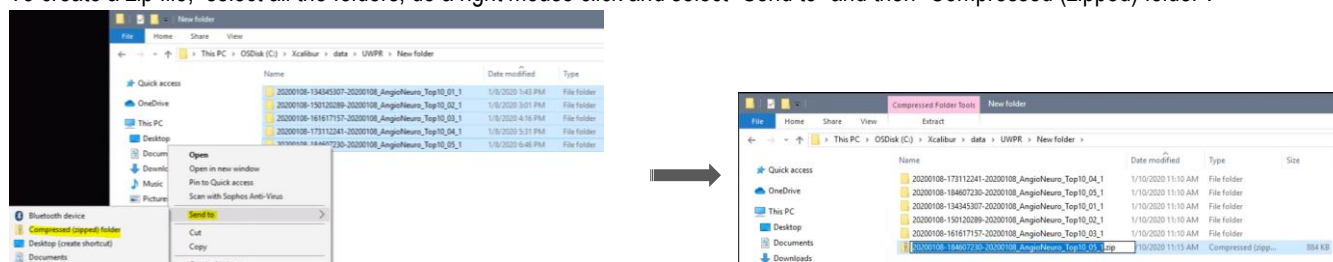


You can rename it or leave the default:

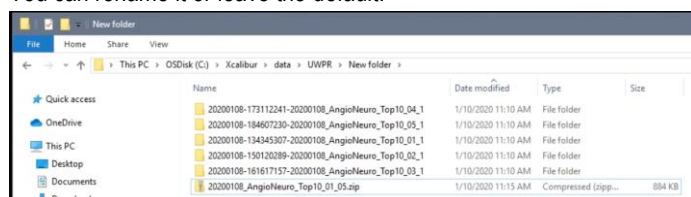


Open .zip file with LogViewer

To open LogViewer (Note this software is supplied by Thermo), double click on the icon on the desktop:

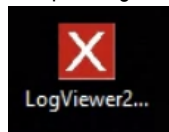


Go to file open, and navigate to the zip file you just created, select it and hit "open":

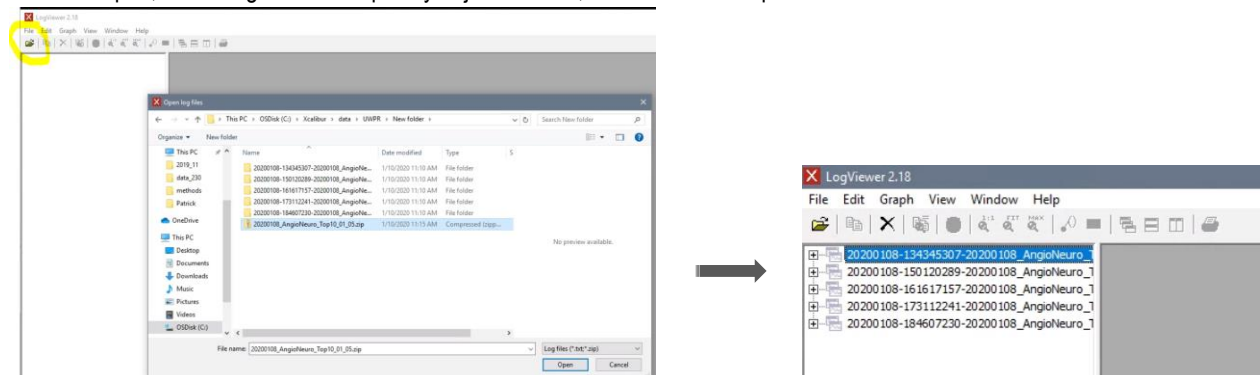


Expand the tree on the left by clicking on the "+", and select the traces you want to look at:



Note it'll likely only show part of the chromatogram, click on "FIT" to see the entire chromatogram:

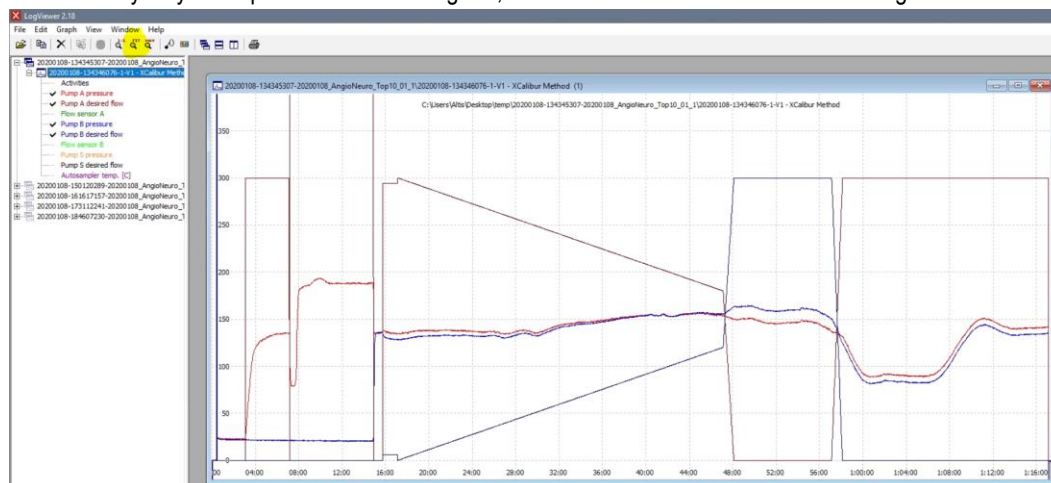

You can open a second one and tile them:

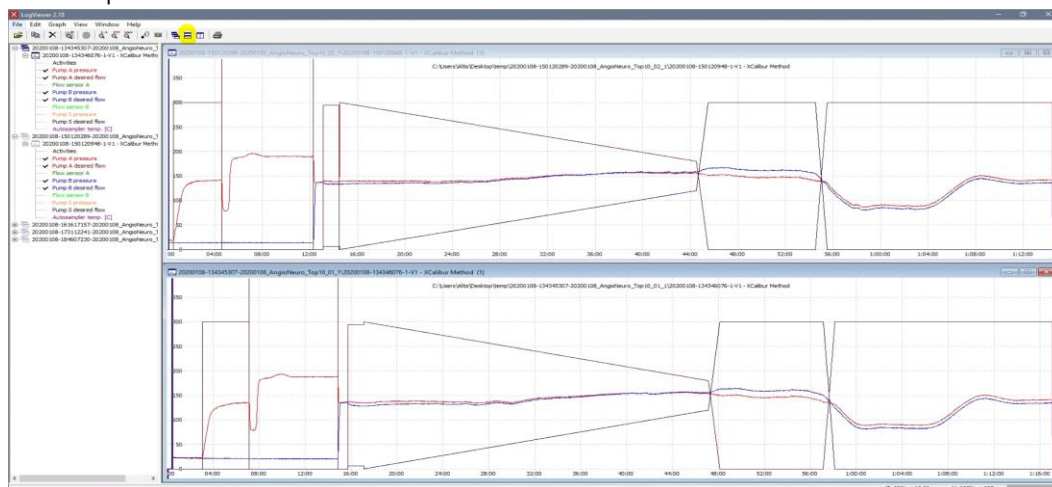


Or you can select “compare graphs” to show them all in one window:



Mount the EASYnLC hard drive to file explorer:

From the start menu, under windows system Open the run command and enter the IP of the EASYnLC in this form \\172.16.0.104\

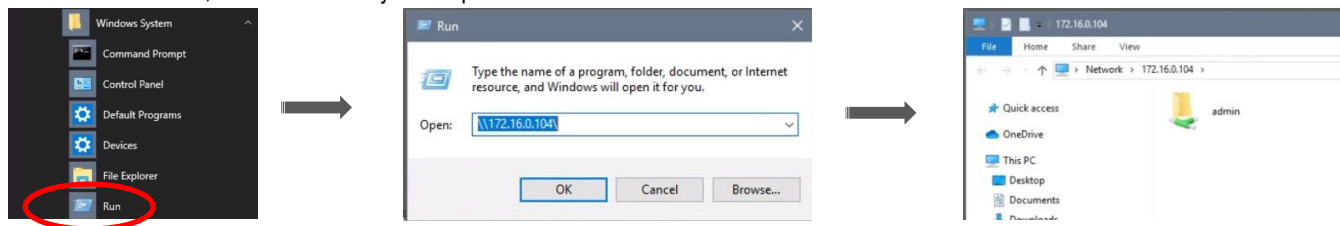

Do a right mouse click on the admin drive and select Map network drive and hit finish:

